Human IL-17A/F ELISpot Kit – Pre-coated

Instructions for use

Catalogue Numbers:

Regular	Pre-coated Plates		
1x96 tests	856.171.001PC		
2x96 tests	856.171.002PC		
5x96 tests	856.171.005PC		

EasySplit	Pre-coated strip plates
1x96 tests	856.171.001PCS
2x96 tests	856.171.002PCS
5x96 tests	856.171.005PCS

For research use only

Fast Track Your Research.....

Table of Contents

1.	Intended use	2
2.	Introduction	2
2.1.	Summary	2
2.2.	Principle of the method	3
3.	Reagents provided	4
4.	Materials/Reagents required but not provided	4
5.	Storage Instructions	4
6.	Safety & Precautions for use	5
7.	Reagent Preparation	5
7.1.	1X Phosphate Buffered Saline (PBS)	5
7.2.	1% BSA PBS Solution (Dilution Buffer)	5
7.3.	Detection Antibody	5
7.4.	Streptavidin – AP conjugate	6
7.5.	BCIP/NBT	6
8.	Sample and Control Preparation	6
8.1.	Cell Stimulation	6
8.2.	Positive Assay Control, IL-17A/F production	6
8.3.	Negative Assay Control	6
8.4.	Sample	7
9.	Method	7
10.	Performance Characteristics	8
10.1	.Specificity	8
10.2	2.Reproducibility and Linearity	8
11.	Bibliography	8

Human IL-17A/F ELISpot Kit – Pre-coated

1. Intended use

Diaclone **ELISpot** is a highly specific immunoassay for the analysis of cytokine and other soluble molecule production and secretion from T-cells at a single cell level in conditions closely comparable to the *in-vivo* environment with minimal cell manipulation. This technique is designed to determine the frequency of cytokine producing cells under a given stimulation and the comparison of such frequency against a specific treatment or pathological state. The ELISpot assay constitutes an ideal tool in the investigation of Th1 / Th2 responses, vaccine development, viral infection monitoring and treatment, cancerology, infectious disease, autoimmune diseases and transplantation.

Utilising sandwich immuno-enzyme technology, Diaclone ELISpot assays can detect both secreted cytokines and single cells that simultaneously produce multiple cytokines. Cell secreted cytokines or soluble molecules are captured by coated antibodies avoiding diffusion in supernatant, protease degradation or binding on soluble membrane receptors. After cell removal, the captured cytokines are revealed by tracer antibodies and appropriate conjugates.

This kit has been configured for research use only and is not to be used in diagnostic procedures.

2. Introduction


2.1. Summary

IL-17A and IL-17F are members of the IL-17 family that play crucial roles in allergic inflammation. They share the strongest homology and the two genes are located on the same chromosome region. Recent studies have shown that activated human CD4+T cells, termed Th17 cells, not only express IL-17A and IL-17F homodimers, but also an IL-17A/F heterodimer cytokine. IL-17A/F is a biologically active protein that induces chemokine production and airway neutrophilia with intermediate potency between IL-17A (most potent) and IL-17F (least potent).


2.2. Principle of the method

A capture antibody highly specific for the analyte of interest is coated to the wells of a PVDF bottomed 96 well microtiter plate either during kit manufacture or in the laboratory. The plate is then blocked to minimise any non-antibody dependent unspecific binding and washed. Cell suspension and stimulant are added and the plate incubated allowing the specific antibodies to bind any analytes produced. Cells are then removed by washing prior to the addition of Biotinylated detection antibodies which bind to the previously captured analyte. Enzyme conjugated streptavidin is then added binding to the detection antibodies. Following incubation and washing substrate is then applied to the wells resulting in coloured spots which can be quantified using appropriate analysis software or manually using a microscope.


1. 96-PVDF bottomed-well plates are first treated with 35% ethanol and then coated with capture antibody.

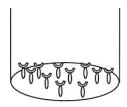

2. Cells are incubated in the presence of the stimulating agent. Upon stimulation they release cytokines which bind to the capture antibodies.

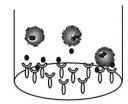


3. Cell removal by washing. Incubation with biotinylated detection antibody.

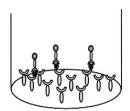

unbound detection Any excess antibodies is removed by washing. Incubation with streptavidin - alkaline phosphatase conjugate.






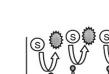


5. Any excess unbound Strep-AP is removed by washing. Incubation with BCIP/NBT.

Finally BCIP/NBT reduction by alkaline phosphatase give a precipitated product which give blue/purple spots.

One spot correspond to one single producing cell.

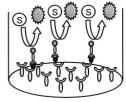


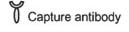






Biotinylated detection antibody

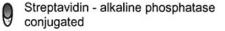



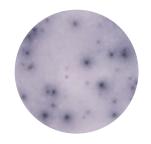

3. Reagents provided

Reagents	.001*	.002*	.005	Reconstitution
Pre-coated 96-well PVDF bottom plates	1	2	5	Rehydrate with 100 µl of PBS1X (see section 9)
Biotinylated Detection antibody	1 (100 µl)	1 (200 µl)	1	Reconstitute with 0.55 ml of distilled water Dilute prior to use (see Detection Antibody, section 7.4)
Streptavidin-Alkaline Phosphatase conjugate	1 (10 µl)	1 (20 µl)	1 (50 µl)	Dilute prior to use (see Streptavidin-AP conjugate, section 7.5)
Bovine Serum Albumin (BSA) – 2 g	1	1	1	Dissolve to prepare dilution buffer (see 1%BSA PBS solution, section 7.3)
Ready to use BCIP/NBT - (Substrate buffer)	1 (11 ml)	1 (25 ml)	2 (25 ml)	Ready to use

^{*}Please note for 001 and 002: detection antibody is provided in liquid form.

4. Materials/Reagents required but not provided

- Miscellaneous laboratory plastic and/or glass, if possible sterile
- Cell culture reagents (e.g. RPMI-1640, L-glutamine, FCS)
- Cell stimulation reagents (e.g. PMA, Ionomycin)
- CO₂ incubator
- Phosphate Buffered Saline (PBS)

5. Storage Instructions

Store kit reagents between 2 and 8°C. Immediately after use remaining reagents should be returned to cold storage (2 to 8°C). Expiry of the kit and reagents is stated on box front labels. The expiry of the kit components can only be guaranteed if the components are stored properly, and if in the case of repeated use of one component, the reagent is not contaminated by the first handling.

6. Safety & Precautions for use

- For research use only not to be used as a diagnostic test.
- Handling of reagents, blood specimens, PBMC, human cell lines should be in accordance with local safety procedures, e.g. CDC/NIH Health manual: "Biosafety in Microbiological and Biomedical Laboratories" 1984.
- Do not eat, drink, smoke or apply cosmetics where kit reagents are used.
- Do not pipette by mouth.
- When not in use, kit components should be stored refrigerated or frozen as indicated on vials or bottles labels.
- All reagents should be warmed to room temperature before use.
- Cover or cap all reagents when not in use.
- Do not mix or interchange reagents between different lots.
- Do not use reagents beyond the expiration date of the kit.
- Use a clean disposable plastic pipette tip for each reagent, standard, or specimen addition in order to avoid cross contamination.
- Use a clean plastic container to prepare the washing solution.
- Thoroughly mix the reagents and samples before use by agitation or swirling.
- All residual washing liquid must be drained from the wells by efficient aspiration or by decantation followed by tapping the plate forcefully on absorbent paper. Never insert absorbent paper directly into the wells.
- When pipetting reagents, maintain a consistent order of addition from well-to-well. This will ensure equal incubation times for all wells.
- **BCIP/NBT substrate** may cause an allergic skin reaction, caution should be taken when handling this reagent, always wear gloves.
- Follow incubation times described in the assay procedure.

7. Reagent Preparation

7.1. 1X Phosphate Buffered Saline (PBS)

For 1 litre of 10X PBS, weigh-out: 80g NaCl

2q KH₂PO₄

14.4g Na₂HPO₄; 2H₂O.

Add distilled water to 1 litre.

Dilute the solution to 1X before use.

Check the pH of the 1X solution and adjust to required pH: 7.4 +/- 0.1.

7.2. 1% BSA PBS Solution (Dilution Buffer)

For one plate, dissolve 0.2 g of BSA in 20 ml of PBS 1X.

7.3. Detection Antibody

Reconstitute the lyophilised antibody with 0.55 ml of distilled water. Gently mix the solution and wait until all the lyophilised material is back into solution.

Please note for 001PC and 002PC kits, detection antibody is provided in liquid form.

If not used within a short period of time, reconstituted Detection Antibody should be aliquoted and stored at -20°C. In these conditions the reagent is stable for at least one year. For optimal performance prepare the reconstituted antibody dilution immediately prior to use.

For one plate, dilute 100 µl of antibody into 10 ml of Dilution Buffer and mix well.

To avoid nonspecific background, it is recommended to filter the working solution using a disposable syringe and a 0.2µm filter disc.

7.4. Streptavidin – AP conjugate

For optimal performance, prepare the Streptavidin-AP dilution immediately prior to use. It is recommended to centrifuge the vial for a few seconds to collect all the volume at the bottom.

For one plate, dilute 10 µl of Streptavidin-AP conjugate into 10 ml of Dilution Buffer and mix well.

Do not keep this solution for further experiments.

To avoid nonspecific background, it is recommended to filter the working solution using a disposable syringe and a 0.2µm filter disc.

7.5. BCIP/NBT

The reagent is ready-to-use.

It should be clear to pale yellow. If precipitates occur, filter the solution using a disposable syringe and a 0.2µm filter disc.

8. Sample and Control Preparation

8.1. Cell Stimulation

Cells can either be stimulated directly in the antibody coated wells (Direct) or, first stimulated in 24 well plates or flask, harvested, and then plated into the coated wells (Indirect).

The method used is dependent on 1) the type of cell assayed 2) the expected cell frequency. When a low number of cytokine producing cells are expected it is also advised to test them with the direct method, however, when this number is particularly high it is better to use the indirect ELISpot method.

All the method steps following stimulation of the cells are the same whatever the method (direct/indirect) chosen.

8.2. Positive Assay Control, IL-17A/F production

We recommend using the following polyclonal activation as a positive control in your assay.

Dilute CD4+ cells in culture medium (e.g. RPMI 1640 supplemented with 2mM L-glutamine and 10% heat inactivated fetal calf serum) containing 1 ng/ml PMA and 500 ng/ml ionomycin (Sigma, Saint Louis, MO). Distribute 1x10⁵ to 2x10⁵ cells per 100 µl in required wells of an antibody coated 96-well PVDF plate and incubate for 15-20 hours in an incubator.

For other stimulators incubation times may vary, depending on the frequency of cytokine producing cells, and should be optimised in each situation.

8.3. Negative Assay Control

Dilute CD4+ cells in culture medium to give an appropriate cell number (same number of unstimulated cells as stimulated sample cells) per 100 µl with no stimulation.

8.4. Sample

Dilute CD4+ cells in culture medium and stimulator of interest (i.e. Sample, Vaccine, Peptide pool or infected cells) to give an appropriate cell number per 100 µl.

Optimal assay performances are observed between 1x10⁵ and 2.5x10⁵ cells per 100 μl.

Stimulators and incubation times can be varied depending on the frequency of cytokine producing cells and therefore should be optimised by the testing laboratory.

9. Method

Prepare all reagents as shown in section 7 and 8.

Note: For optimal performance prepare the Streptavidin-AP dilution immediately prior to use.

Assay Step	Details				
Addition	Add 100 μl of PBS 1X to every well				
Incubation	Incubate plate at room temperature (RT) for 10 min				
Wash	Empty the wells by flicking the plate over a sink & gently tapping on absorbent paper				
Addition	Add 100 µl of sample , positive and negative controls cell suspension to appropriate wells providing the required concentration of cells and stimulant (cells may have been previously stimulated see section 8.)				
Incubation	Cover the plate and incubate at 37°C in a CO ₂ incubator for an appropriate length of time (15-20 hours) Note: do not agitate or move the plate during this incubation				
Addition	Empty the wells and remove excess solution then add 200 µl of PBS 1X to every well				
Incubation	Incubate the plate at 4°C for 10 min				
Wash	Empty the wells as previous and wash the plate 3x with 200 µl of PBS 1X				
Addition	Add 100 μl of diluted detection antibody to every well				
Incubation	Cover the plate and incubate at RT for 1 hour 30 min				
Wash	Empty the wells as previous and wash the plate 3x with 200 μl of PBS 1X				
Addition	Add 100 µl of diluted Streptavidin-AP conjugate to every well				
Incubation	Cover the plate and incubate at RT for 1 hour				
Wash	Empty the wells and wash the plate 3x with 200 μl of PBS 1X				
Wash	Peel off the plate bottom and wash both sides of the membrane 3x under running distilled water, once washing complete remove any excess solution by repeated tapping on absorbent paper				
Addition	Add 100 μl of ready-to-use BCIP/NBT buffer to every well				
Development	Incubate the plate for 5-15 min monitoring spot formation visually throughout the incubation period to assess sufficient colour development				
Wash	Empty the wells and rinse both sides of the membrane 3x under running distilled water. Completely remove any excess solution by gentle repeated tapping on absorbent paper the wells to dry and then read results. The frequency of the resulting coloured spots				
	Incubation Wash Addition Incubation Addition Incubation Wash Wash Wash Wash Addition Development Wash				

Read Spots: allow the wells to dry and then read results. The frequency of the resulting coloured spots corresponding to the cytokine producing cells can be determined using an appropriate ELISpot reader and analysis software or manually using a microscope.

Note: spots may become sharper after overnight incubation at 4°C in the dark

Plate should be stored at RT away from direct light, but please note that colour may fade over prolonged periods so read results within 24 hours.

10. Performance Characteristics

10.1. Specificity

The assay recognizes natural human IL-17A/F.

To define specificity of this IL-17A/F antibody pair, several proteins were tested for cross reactivity. There was no cross reactivity observed for any of these proteins tested : IL-17A, IL-17B, IL-17D, IL-17E, IFN γ , Gal-1 and IL-33). This antibody pair shows cross reactivity with human IL-17F. This testing was performed using the equivalent human IL-17A/F antibody pair in an ELISA assay.

10.2. Reproducibility and Linearity

Intra-assay reproducibility and linearity were evaluated by measuring the spot development following the stimulation (PMA / Ionomycin) of 6 different CD4+ cell concentrations, 12 repetitions. The data show the mean spot number, range and CV for the six cell concentrations.

Cells / well	n	Mean number of spots per well	Min	Max	CV%
200000 recommended	12	43	33	55	16
100000 recommended	12	47	40	56	11
50000	12	33	21	65	36
25000	12	17	12	26	28
12500	12	15	8	23	28
6250	12	11	4	18	41

11. Bibliography

- 1. Chang S.H. et.al: A novel heterodimeric cytokine consisting of IL-17 and IL-17F regulates inflammatory responses. Cell Research (2007) 17:435-440.
- 2. Kawaguchi et.al: IL-17 cytokine family. J Allergy Clin Immunol (2004) vol 114, n 6.
- 3. Wright et.al: Identification of an interleukin 17F/17A Heterodimer in Activated Human CD4+ T cells. Journal of Biological Chemistry (2007) vol.282, NO.18, pp. 134447-13455.
- 4. Wright et.al: The Human IL-17F/IL-17A Heterodimeric Cytokine Signals through the IL-17RA/IL-17RC Receptor Complex. J. Immuno, 2008, 181: 2799-2805.
- 5. Cheung et.al: Molecular Mechanisms of Cytokine and Chemokine Release from Eosinophils Activated by IL-17A, IL-17F, and IL-23: Implication for Th17 Lymphocytes-Mediated Allergic Inflammation. J.Immuno. 2008. 180:5625-5635.

Products Manufactured and Distributed by:

Diaclone SAS 6 Rue Dr Jean-François-Xavier Girod 25000 Besançon France Tel +33 (0)3 81 41 38 38

Email: techsupport@medixbiochemica.com